



# BVHDWRSMSGGA

## Bella Vita Handheld Wireless Massager


### Quick Start Guide

## Safety & Warnings

- Read all instructions before using the appliance and retain for reference.
- Disconnect the power adaptor before cleaning.
- Do not use water to clean.
- Cease usage immediately if there are any irregularities in operation, including electrical fault or burning smell, and contact Kogan Support.
- Do not remove the power cord with a wet hand.
- Do not use in humid areas, such as the bathroom.
- Do not use any coating, volatile oil or pesticide with the appliance.
- Unplug when not in use.
- Do not allow the product to be exposed to sunlight for a long time.
- This appliance is not a toy. Do not use near children.

## Overview

- Insert the plug into the socket and press the power switch: 5 modes and 5 speed are selectable.
- Massaging 1-2 times daily for 5-15 minutes can be therapeutic.
- The infrared heat function raises blood circulation.
- The appliance will switch off automatically after 20 minutes.

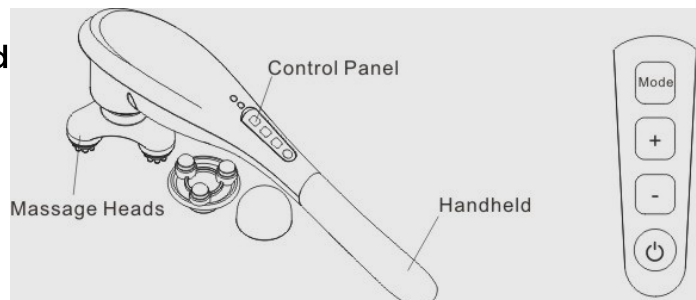

## Specifications

Rated Voltage: DC9V  
Rated Wattage: 9W  
Vibration Partial Rate: 5mm  
Vibration Frequency: approx. 4000 RPM

## Cleaning & Care

- Use a neutral cleaner and soft cloth to wipe.
- Do not use harsh cleaners or allow any liquid to enter the device.

User manual is subject to change without notice. For the latest version of your user manual, please visit <https://www.kogan.com/usermanuals/>

## Attention

Please handle this product with care and inspect it regularly to ensure it is in good working order.

If the product, power supply cord or plug shows any signs of damage: stop use, unplug and contact [Kogan.com](http://Kogan.com) support.