

# **ELECTRIC TOOTHBRUSH WITH UV SANITISER KATBRSHVSA QUICK START GUIDE**

## **IMPORTANT SAFEGUARDS**

This electric toothbrush is safe and effective when used in accordance with the directions. Like all electrical appliances however, precautions must be followed.

### **Danger - Water**

To reduce the risk of burns, electrocution, fire or injury to persons:

- Do not place or store the charger base with UV sanitiser where it can fall or be pulled into a bathtub, sink or toilet.
- Do not place or drop the charger into water or any other liquids.
- Do not reach for a product that has fallen into water or another liquid. Unplug it from the power socket if it is safe to do so.
- Do not use this product while bathing. Unplug the charger base before cleaning, and ensure that it is dry before plugging it back into an electrical outlet.

### **Danger - Voltage**

- Do not plug the charger base into an outlet with a voltage other than the one specified on the bottom of the charger.
- This product is designed to operate within a range of 100-240 volts.
- Voltage converters do not guarantee compatibility.
- Do not force the plug into an outlet.
- If the plug does not fit easily into the outlet, do not use it.
- Never wrap the cord tightly around the appliance. Never allow the power supply cord to be pulled, twisted or severely bent. Untwist the cord before use.

### **Danger - Warnings**

- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given adequate supervision or instruction concerning use of the product by an adult responsible for their safety.
- Children should be supervised to ensure they do not play with the product.
- Use this product for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only.
- Do not use any attachments or brush heads other than those recommended by the manufacturer.
- The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution. If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid damaging the plastic.
- Never operate this product if it has a damaged cord, plug or charger with UV sanitiser.
- Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- Keep the product, including power cord away from heat surfaces.
- If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- Do not attempt to remove the battery
- Do not clean the toothbrush in the dishwasher.


- Discontinue use of the charger with UV sanitiser if the UV light bulb remains on when the lid is open.
- UV light can be harmful to human eyes and skin. Keep the charger with UV sanitiser out of the reach of children at all times.
- The UV light bulb is hot during or after the sanitising cycle. Do not touch the UV light bulb when it is hot.
- Do not look directly into the UV light.
- Do not allow any impact to the UV lamp. If the UV lamp is cracked, chipped or damaged, do not operate.
- In the event of the UV lamp breaking, leave the location for 20-30 minutes to allow the fumes to settle, and then clean it with thick paper towelling. Do not touch the mercury directly with bare hands.
- For the best disinfection performance, always keep the UV lamp clean and avoid stains.
- In case the UV lamp gets dirty, wipe it with a soft cloth very gently.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Unplug from outlet when you will not be using the charger for a long period of time, and before servicing.
- Do not pull or carry by cord, use cord as a handle, close a door on cord, or pull cord around sharp edges or corners.
- Do not run appliance over cord.
- Do not handle plug or appliance with wet hands.
- Turn off all controls before unplugging.
- Replacement of the UV lamp can only be carried out by a Kogan technician. It should not be replaced by the user.
- The UV cover automatically switches off after the sanitiser cycle to reduce the risk of excessive ultraviolet radiation.

## MEDICAL WARNINGS


- If you have had oral surgery or gum surgery within the past two months, consult your dental provider before using the product.
- Consult your dental provider if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
- Consult your physician prior to using the toothbrush if you have any medical concerns.
- Never use while sleeping or drowsy.

## OVERVIEW OF THE PRODUCT


Handle

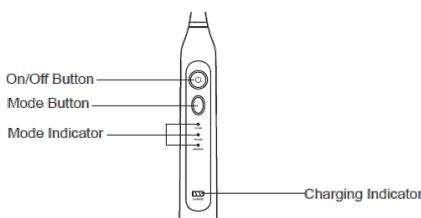

Charger

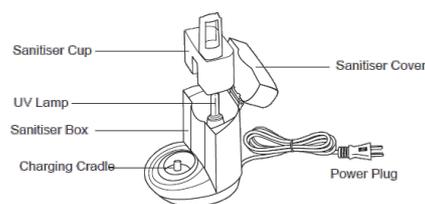

Brush Heads

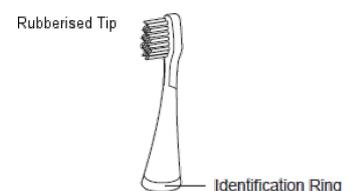


# GETTING STARTED

Charging the toothbrush: Be sure to turn the toothbrush off before putting it into the charger. Plug the power plug into a wall outlet, then place the toothbrush handle into the charger. The charger indicator will light up, indicating the handle is charging. Once the battery is fully charged, all 3 sections of the charging indicator will be lit up in solid green, indicating the toothbrush is fully charged.


Note:

- It takes at least 24 hours to charge the battery fully. When fully charged, the toothbrush can be used for 20 pre-timed brush cycles (approximately 10 days) if brushing twice daily.
- The toothbrush must be charged for 24 hours if:
  - this is the first time using the toothbrush
  - the toothbrush has not been used for more than 3 months.
  - the charge is fully drained from use
- To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.
- During brushing, the toothbrush will get warm. This is completely normal.
- The operating duration will become shorter as the battery gets older. If the toothbrush operates for only 1 short time after being charged, then this means that the battery is nearing the end of its life cycle. Contact the Kogan support center for assistance.
- Attach a brush head to the handle. There will be about 0.5mm space between the brush head and the handle. (refer to Fig 1)
- Press the mode button to select the mode. If you know the mode used last time, then there is no need to select again. There are 3 modes:
  - Clean: standard mode for superior teeth polishing
  - Polish: Teeth cleaning with polishing function
  - Massage: Pulsates for gentle gum stimulation.
- Wet the bristles and apply a small amount of toothpaste.
- Place the bristles against the teeth at a slight angle towards the gum line. (refer to Fig 2)
- Press the on/off button to turn the toothbrush on.
- In order to maximise the cleaning effectiveness, divide your mouth into 4 sections. Start with section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when it is time to move to the next section. The toothbrush will automatically stop after 2 minutes use. (refer to Fig 3)
- Take the brush head away from the handle. (refer to Fig 4)
- Rinse both the brush head and handle under running water and then wipe them dry. (refer to Fig 5)
- Put the handle back into the charger.

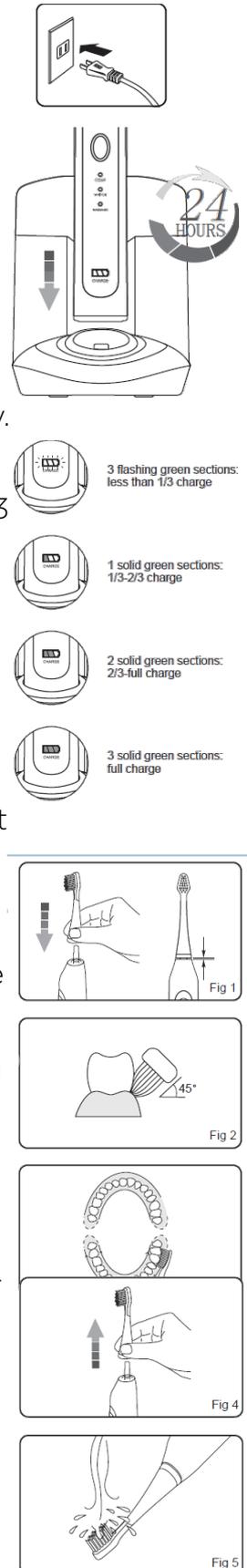

## SANITISING THE BRUSH HEAD

This product includes a UV sanitiser to sanitise the brush heads between use. The sanitiser has proven to be 99.9% effective against E.coli, Candida Albicans, Salmonella and H. Simplex.

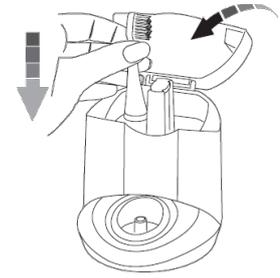


1. Gently shake the toothbrush head to remove any excess water after cleaning.
2. Open the sanitiser cover.
3. Place the brush heads on the stand with the bristles directly facing the UV lamp.
4. Make sure the charger with the UV sanitiser is plugged into a working outlet.
5. Close the sanitiser cover to activate the UV sanitising process.
  - (a) The sanitiser cover will glow with a slight blue hue indicating the UV lamp is lit and the sanitising cycle has started. Please note that this glow may be hard to see in well-lit areas.
  - (b) When the sanitising cycle is complete, the sanitiser will automatically turn off and the cover will no longer glow.
6. When using the sanitiser:
  - (a) Place only compatible brush heads in the sanitiser.
  - (b) The UV sanitiser will run for approximately seven minutes, then shut off automatically.
  - (c) The UV sanitiser will run automatically every time the UV sanitiser cover is closed.
  - (d) Do not place other objects, including other model brush heads, into the UV sanitiser.

## CLEANING THE UV SANITISER

Be sure to remove the plug from the electrical outlet before cleaning the product.

1. Open the sanitiser cover to remove.
2. Using a soft cloth, wipe off any dirt on the charger and UV lamp.
3. Clean exterior surfaces and plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheating or fire.
4. Do not clean the sanitiser, charger or brush heads in the dishwasher.

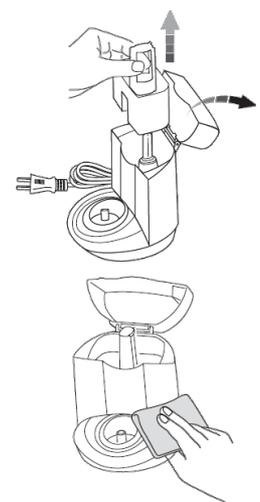

# TROUBLESHOOTING

<b>Problem</b>	<b>Possible Cause</b>	<b>Solution</b>
Toothbrush doesn't work	Product has just been purchased, or has not been used for more than 3 months	Charge the handle for at least 24 hours
Weak or diminished vibrations	Toothbrush is pressed too hard against your teeth	Place brush head gently against your teeth.
Toothbrush operates for a short period only, even after being charged	Battery has reached the end of its life cycle	Contact the Kogan customer support team
	Charging time was insufficient	Charge the handle for at least 24 hours.
UV sanitiser light does not turn on	Sanitiser cover may be open	Close the cover
	Power plug may not be plugged in	Plug the power plug in
	UV lamp may be defective	Contact the Kogan customer support team